



Hillrom™

## Explore Welch Allyn Green Series 777 Integrated Wall Systems & Transformers featuring our latest Welch Allyn Exam Tools

Welch Allyn **Green Series 777** Integrated Wall Systems and Transformers can help you stay efficient during your examinations. Choose our wide-view **PanOptic Plus** Ophthalmoscope and **MacroView Plus** Otoscope, digital or manual blood pressure devices and integrated or stand-alone oral/axillary/rectal thermometry with a Welch Allyn **SureTemp Plus** Thermometer.

Choose wide-view, digitally-capable devices or standard heads.

Turns on automatically when handle is removed from cradle and turns off when you return it.

Minimize cross-contamination risk and dispense a new disposable ear tip with each exam.

Choose a stand-alone BP device and Thermometer or order a compatible Welch Allyn multi-parameter vital signs monitor.




- Easy to upgrade devices from the Welch Allyn 767-series to the 777 Integrated Wall System
- Automatic on/off sensors help conserve energy
- Handle-based rheostats for quick light intensity adjustments
- Modular coiled cords designed to withstand daily use

## SPECIFICATIONS

<b>Shipping Dimensions &amp; Weight</b>	Width: 34.33 in; Height: 12.83 in; Depth: 9.48 in; Weight: 24.25 lbs
<b>Installed Dimensions</b>	Width: 30.00 in; Height: 11.75 in
<b>Coiled Power Cord Length (extended)</b>	12 ft
<b>Operating Environment</b>	Temperature Limitation: 15°C to 40°C (59°F to 104°F) Humidity Limitation: 15% to 95% Atmospheric Limitation: 700 hPa to 1060 hPa
<b>Transport &amp; Storage</b>	Temperature Limitation: -20°C to 50°C (-4°F to 120°F) Humidity Limitation: 15% to 95% Atmospheric Limitation: 700 hPa to 1060 hPa
<b>Safety Class</b>	II Intermittent Operation
<b>Electrical</b>	Input: 100–240 VAC, 50–60 Hz, 0.18A Output: 5 VDC 1A
<b>Materials</b>	ABS and elastomer rubber. Made from non-latex material.
<b>Approvals</b>	IEC 60601-1, UL 60601-1, IEC 60601-1-2, CAN/CSA C22.2 No 601.1-M90, AS/NZS3200

## GUIDANCE ON INSTALLATION

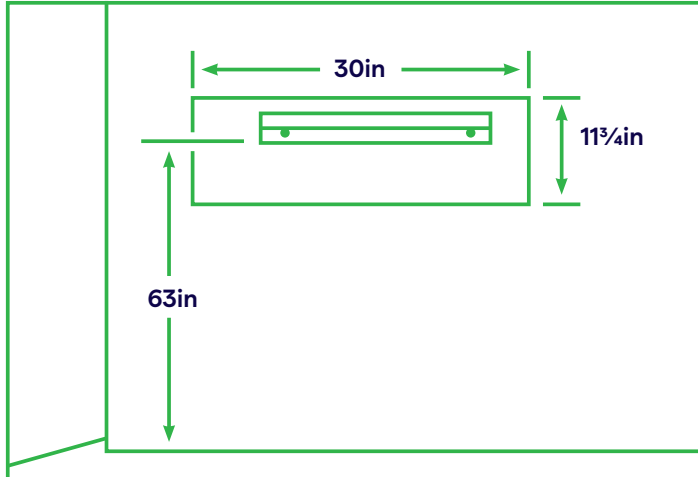

## PARTS & ACCESSORIES

- **118-3-US: PanOptic Plus Ophthalmoscope** (head only)
- **238-3: MacroView Plus Otoscope** (head only)
- **21504: Insufflator Bulb for MacroView Otoscope**
- **11840-IEX-US: iExaminer SmartBracket Accessory for MacroView Plus Otoscope**
- **52432-CLR-1: 2.75 mm pediatric size LumiView Clear Single-Use Ear Specula**
- **52434-CLR-1: 4.25 mm adult size LumiView Clear Single-Use Ear Specula**
- **05031-750: SureTemp Plus Thermometer Disposable Probe Covers (Case of 7,500)**
- **104299: Coiled Cord and Handle Assembly (12 ft)**
- **PWCD-GS777-B: Green Series 777 Wall System & Transformer Power Cord Kit**
- **772320: Green Series 777 Wall System & Transformer USB Power Cord (18 in)**

Connect with your **Baxter** representative today to discuss the best **Welch Allyn Physical Exam & Diagnostic** solution for your clinical needs.

**hillrom.com**

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

Baxter, Connex, Green Series, Hillrom, iExaminer, LumiView, MacroView, PanOptic, SmartBracket, SureTemp, and Welch Allyn are trademarks of Baxter International Inc. or its subsidiaries. ALL RIGHTS RESERVED. APR348801 Rev 1 27-OCT-2022 ENG-US